



HORIZON 4.0 PRECISION PAEDIATRIC INTENSIVE CARE BED

**AT LAST A PICU BED WITH 360° ACCESS AND A
HIGH PERFORMANCE INTEGRATED WEIGH
SCALE ACCURATE TO 50 GRAMS AND FULL
LENGTH X - RAY CAPABILITY FROM
UNDERNEATH**

Favero Health Projects

The Horizon 4.0 Precision PICU bed features an integrated approved weigh scale. There is a tare function and an add/take away function. The screen is simple and intuitive to use. It has a function to set minimum and maximum thresholds in order to check patient activity and is accurate to 50 grams all the way up to 150 kg. This allows accurate administration of medications.



In intensive care you have told us that very important to have access from the head end in case a patient needs to be intubated, the removable headboard panels are designed with this in mind. Furthermore, lowering the side rails allows 360° access to the child, which is essential to prevent incidents with ventilators or cables and therefore guarantees maximum safety to the patient at all times.



Address: 31 Airlie Road, Pullenvale, Queensland 4069 Australia
Phone: 07 3166 9769 Fax: 07 3011 1029 Mobile: 0410 537 746
Email: sales@newmedical.com.au Web: www.newmedical.com.au
ABN: 63 124 888 348



Radio-translucent sleep deck sections in HPL. More resistant and durable. X-Ray compatible.



Potter guides compatible with both traditional and digital x-ray cassettes. X Ray along almost the entire length of the underside.



NEW COMMAND SYSTEM





10 Key Operator Handset



Operator Control Panel

NEW COMMAND SYSTEM

- Visual command activation display.
- Back rest and sleep deck angle indicators.
- Battery and command lock indicator display.
- One touch emergency buttons.
- Synchronized one touch bed commands via the +FUNC key including AUTOCONTOUR, EXAM, FOWLER, SEAT AND FAINT positions.
- Harmonised and simultaneous command interaction between handset and control panel.

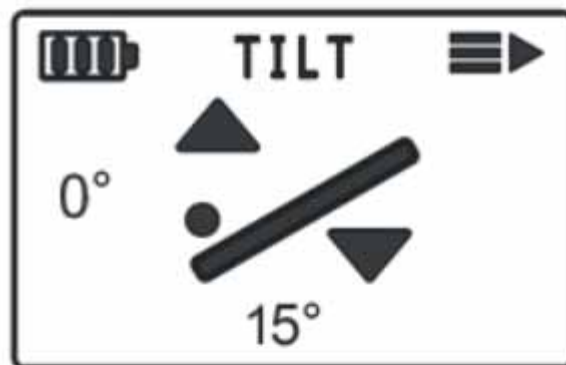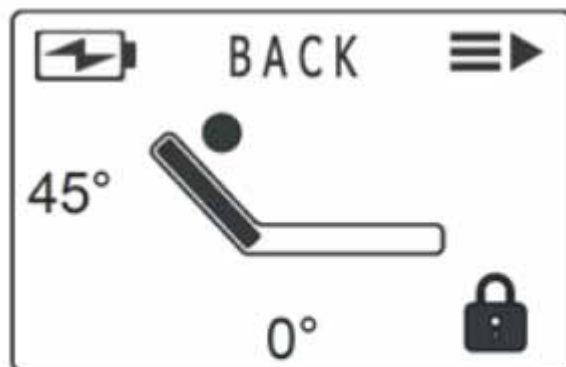

TECHNICAL DATA SHEET

Technical Data System	
Platform	4 Sections
Mattress Size	1314 x 655 mm
Platform Size	1314 x 655 mm
Footprint	1574 x 930 mm
Height of Mattress base from floor	From min. 640 to max 870 mm
Location of the Operator Control Panel	At foot end and additional at head end
Height of Siderails from platform	378 mm
Safe Working Load	90 kg
Wheels	4 150 mm Twin Wheels with Braking System
Backrest Inclination	0° to 75° Gas Spring Assisted
Backrest Inclination	0° to 30°
Trendelenburg position	Max 15°
Reverse Trendelenburg position	Max 15°
Weighting System	Yes

Accessories	
9LP9201	Paediatric critical care mattress 1314 x 655 x 100 mm, with waterproof, breathable, zippable cover.
9LP9036	Pull-out standard accessories bar at the head end 25 x 10mm
9LP9037	ECMO Support

